




**CELLular**  
**Dynamics**  
international

# **iCell<sup>®</sup> Astrocytes Prototype User's Guide**



## Trademarks

iCell is a registered trademark, and Cellular Dynamics and the  logo are trademarks of Cellular Dynamics International, Inc. and may not be used without the express written permission of Cellular Dynamics International, Inc. (CDI).

All other brands, product names, company names, trademarks, and service marks are the properties of their respective owners.

## Restrictions and Liabilities

This document is provided “as is.” CDI assumes no responsibility for any typographical, technical, or other inaccuracies in this document. CDI reserves the right to periodically change information that is contained in this document; however, CDI makes no commitment to provide any such changes, updates, enhancements, or other additions to this document to you in a timely manner or at all.

THE iCELL ASTROCYTES ARE PROVIDED “AS IS” TO YOU. TO THE FULLEST EXTENT PERMITTED BY APPLICABLE LAW, CDI MAKES NO, AND DISCLAIMS ALL, REPRESENTATIONS, WARRANTIES, CONDITIONS OR COVENANTS, EXPRESS OR IMPLIED, WITH RESPECT TO ANY PRODUCT REFERENCED HEREIN OR PERFORMANCE OF ANY SERVICES REFERENCED HEREIN, INCLUDING BUT NOT LIMITED TO, ANY EXPRESS OR IMPLIED WARRANTIES OR CONDITIONS OF FITNESS FOR A PARTICULAR PURPOSE, NON-INFRINGEMENT, MERCHANTABILITY, DURABILITY, TITLE, OR RELATED TO THE PERFORMANCE OR NON-PERFORMANCE OF ANY PRODUCT REFERENCED HEREIN OR PERFORMANCE OF ANY SERVICES REFERENCED HEREIN).

This document might contain references to third-party sources of information, hardware or software, products or services and/or third-party websites (collectively the “Third-Party Information”). CDI does not control, and is not responsible for, any Third-Party Information, including, without limitation the content, accuracy, copyright compliance, compatibility, performance, trustworthiness, legality, decency, links, or any other aspect of Third-Party Information. The inclusion of Third-Party Information in this document does not imply endorsement by CDI of the Third-Party Information or the third party in any way.

CDI does not in any way guarantee or represent that you will obtain satisfactory results from using iCell Astrocytes as described herein. You assume all risk in connection with your use of iCell Astrocytes.

## Conditions of Use

iCell Astrocytes are for life science research use only and subject to the use restrictions as contained in Appendix A. You are responsible for understanding and performing the protocols described within. CDI does not guarantee any results you may achieve. These protocols are provided as CDI’s recommendations based on its use and experience with iCell Astrocytes.

## Origin

iCell Astrocytes are manufactured in the United States of America.

## Copyright Notice

© 2014 Cellular Dynamics International, Inc. All rights reserved. This document may not be reproduced, distributed, modified or publicly displayed without the express written permission of Cellular Dynamics International, Inc.

## Revision History

Version 0.3: October 2014

Version 0.2: August 2014

Version 0.1: March 2013

# Table of Contents

|  |           |
|--|-----------|
| <b>Before You Begin</b> .....  | <b>ii</b> |
| <b>Chapter 1. Introduction</b> .....   | <b>1</b>  |
| Components Supplied by Cellular Dynamics .....                                     | 2         |
| Required Equipment and Consumables .....   | 2         |
| Technical Support and Training.....  | 3         |
| Workflow Diagram .....   | 4         |
| <b>Chapter 2. Handling and Storage</b> .....                                       | <b>5</b>  |
| <b>Chapter 3. Preparing Cell Culture Surfaces</b> .....                            | <b>6</b>  |
| <b>Chapter 4. Preparing the Medium</b> .....                                       | <b>7</b>  |
| <b>Chapter 5. Thawing iCell Astrocytes</b> .....                                   | <b>8</b>  |
| <b>Chapter 6. Plating iCell Astrocytes</b> .....                                   | <b>10</b> |
| <b>Chapter 7. Maintaining iCell Astrocytes</b> .....                               | <b>11</b> |
| <b>Appendices</b> .....  | <b>12</b> |
| Appendix A. Intellectual Property Rights, Use Restrictions, and Limited License... | 12        |
| Appendix B. Product Provided “AS IS” .....   | 12        |
| Appendix C. Limited Liability.....   | 13        |

## Before You Begin

- Immediately transfer the frozen vials to liquid nitrogen storage.
- Read this entire iCell® Astrocytes Prototype User's Guide before handling or using iCell Astrocytes.
- iCell Astrocytes are for life science research use only. See Appendix A for more information and other restrictions.
- A Safety Data Sheet (SDS) for dimethyl sulfoxide (DMSO), in which iCell Astrocytes are frozen, is available online at [www.cellulardynamics.com/lit/](http://www.cellulardynamics.com/lit/) or on request from Cellular Dynamics International. Only technically qualified individuals experienced in handling DMSO and human biological materials should access, use, or handle iCell Astrocytes.

Notes

## Chapter 1. Introduction

Cellular Dynamics International's (CDI) iCell Astrocytes are a highly pure population of human astrocytes derived from induced pluripotent stem (iPS) cells using CDI's proprietary differentiation and purification protocols. iCell Astrocytes exhibit expected physiological characteristics and responses. These cells provide a reliable source of human astrocytes suitable for use in targeted drug discovery, toxicity testing, and other life science research.

## Components Supplied by Cellular Dynamics

Notes

| Item  | Catalog Number     |
|---|--------------------|
| iCell Astrocytes Prototype <sup>1</sup>   | ASC-100-020-001-PT |
| iCell Astrocytes Prototype User's Guide <sup>1</sup>  |                    |
| Certificate of Testing <sup>2</sup>   |                    |
| Certificate of Origin<br>If required for shipping purposes  |                    |
| <p>1 Safety Data Sheet and User's Guide available online at <a href="http://www.cellulardynamics.com/lit/">www.cellulardynamics.com/lit/</a></p> <p>2 Available by emailing <a href="mailto:support@cellulardynamics.com">support@cellulardynamics.com</a> or calling (877) 320-6688 (US toll-free) or (608) 310-5100</p> |                    |

## Required Equipment and Consumables

| Item   | Vendor            | Catalog Number |
|--|-------------------|----------------|
| <b>Equipment</b>   |                   |                |
| 37°C Water Bath  | Multiple Vendors  |                |
| Biological Safety Cabinet with UV Lamp                   | Multiple Vendors  |                |
| Cell Culture Incubator                                   | Multiple Vendors  |                |
| Hemocytometer or Automated Cell Counter                  | Multiple Vendors  |                |
| Liquid Nitrogen Storage Unit                             | Multiple Vendors  |                |
| Pipettors  | Multiple Vendors  |                |
| Tabletop Centrifuge                                      | Multiple Vendors  |                |
| <b>Consumables</b>                                       |                   |                |
| 15 ml and 50 ml Centrifuge Tubes                         | Multiple Vendors  |                |
| 6-well Cell Culture Plates                               | Nunc              | 140675         |
| Astrocyte Culture Medium*                                | Multiple Vendors  |                |
| DMEM   | Life Technologies | 10569          |
| Growth Factor Reduced Corning Matrigel Matrix (Matrigel) | Corning           | 354230         |
| Pipettes   | Multiple Vendors  |                |

\* Recommended astrocyte culture medium: DMEM, High Glucose, GlutaMAX, Pyruvate (Life Technologies, Cat. No. 10569) with 1X N-2 Supplement (Life Technologies, Cat. No. 17502048) and 10% FBS (Hyclone, Cat. No. SH30396.03) or similar commercially available astrocyte culture medium.

---

## Technical Support and Training

CDI's Technical Support Scientists have the necessary laboratory and analytical experience to respond to your inquiries. In addition, in-lab training may be available upon request.


**Telephone** (877) 320-6688 (US toll-free) / (608) 310-5100 x5  
Monday - Friday, 8:30 am - 5:00 pm US Central Time

**Fax** (608) 310-5101

**Email** [support@cellulardynamics.com](mailto:support@cellulardynamics.com)

## Workflow Diagram

Notes

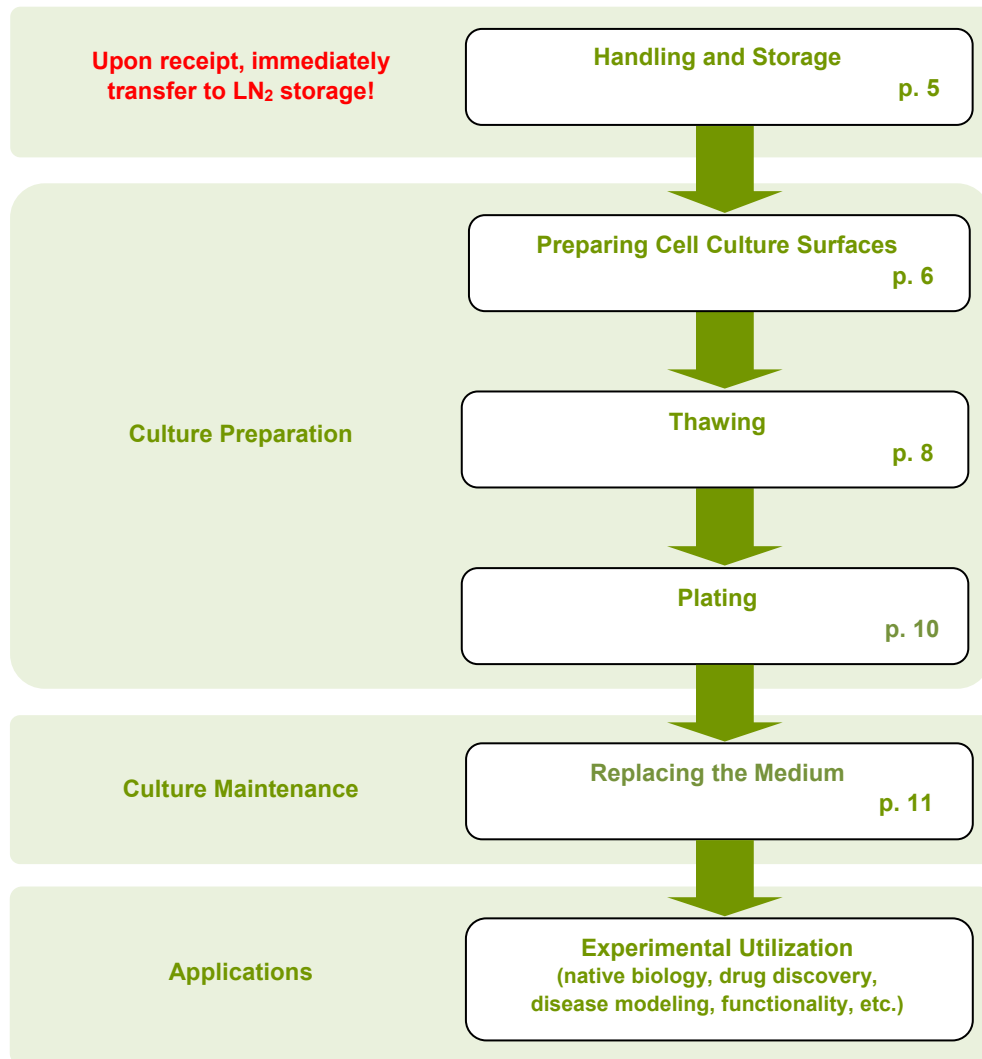



## Chapter 2. Handling and Storage

iCell Astrocytes are provided as cryopreserved single-cell suspensions in 1.5 ml cryovials. Upon receipt, directly transfer the cryobox containing iCell Astrocytes to the vapor phase of a liquid nitrogen storage dewar. CDI strongly recommends transferring the entire cryobox into the storage rack to avoid transferring individual vials.



*It is critical to maintain cryopreserved iCell Astrocytes at a stable temperature. Minimize exposure of cryopreserved iCell Astrocytes to ambient temperature when transferring vials to liquid nitrogen storage.*

## Chapter 3. Preparing Cell Culture Surfaces

iCell Astrocytes will plate and function on cell culture plates pre-coated with Matrigel. The following procedure details coating 6-well cell culture plates. Scale volumes appropriately for other well formats.

1. Dilute Matrigel in ice-cold DMEM to a final concentration of 0.083 mg/ml.
2. Immediately add 1 ml/well of a 6-well cell culture plate.
3. Incubate the plate at room temperature for at least 1 hour before plating iCell Astrocytes.

**Note:** *Plates coated with Matrigel can be stored at 4°C for up to 1 week. Equilibrate the plates in a 37°C cell culture incubator before use.*

Notes

## Chapter 4. Preparing the Medium

CDI recommends thawing, plating, and maintaining iCell Astrocytes in the following medium or in similar commercially available astrocyte culture medium.

| Component                              | Final Concentration |
|--|---------------------|
| DMEM, High Glucose, GlutaMAX, Pyruvate | 89%                 |
| FBS                                    | 10%                 |
| N-2 Supplement (100X)                  | 1X                  |

**Table 1: Components of Astrocyte Culture Medium**

## Chapter 5. Thawing iCell Astrocytes

Maintain iCell Astrocytes in liquid nitrogen until immediately before thawing to ensure maximal performance of the cells. Completing the following steps in a time-efficient manner facilitates optimal iCell Astrocytes viability and performance.

**Note:** Thaw no more than 3 vials of iCell Astrocytes at one time.

1. Equilibrate the astrocyte culture medium at room temperature for 2 - 4 hours before thawing iCell Astrocytes.
2. Remove the iCell Astrocytes cryovial from the liquid nitrogen storage tank.

**Note:** If necessary, place cryovials on dry ice for up to 10 minutes before thawing.

3. Immerse the cryovial in a 37°C water bath for 3 minutes (avoid submerging the cap) holding the tube stationary (no swirling). Use of a floating microcentrifuge tube rack is recommended.
4. Immediately remove the cryovial from the water bath following the 3 minute incubation, spray with 70% ethanol, and place in a biological safety cabinet.
5. Gently transfer iCell Astrocytes cryovial contents to a sterile 50 ml centrifuge tube using a 1 ml pipettor.

**Note:** Use of a 50 ml centrifuge tube facilitates suitable mixing to minimize osmotic shock and increase cell viability.

6. Rinse the empty iCell Astrocyte cryovial with 1 ml of room temperature astrocyte culture medium to recover any residual cells from the cryovial. Transfer the 1 ml astrocyte culture medium rinse from the cryovial drop-wise (~1 drop/second) to the 50 ml centrifuge tube containing the iCell Astrocytes suspension. Gently swirl the tube while adding the medium to mix the solution completely and minimize the osmotic shock on the thawed cells.



Drop-wise addition of astrocyte culture medium to the cell suspension is critical to minimize osmotic shock and ensure maximum viability and subsequent attachment of the cells to the plating substrate.

7. Slowly add 8 ml of room temperature astrocyte culture medium to the 50 ml centrifuge tube drop-wise (~1 - 2 drops/second). Gently swirl the centrifuge tube while adding the medium.



It is critical to add the 8 ml of astrocytes culture medium slowly to ensure maximum viability and attachment of the cells once plated.

## Notes

8. Gently mix the contents of the 50 ml centrifuge tube by swirling or inverting 2 - 3 times. Gentle mixing is critical to ensure maximum viability. Avoid vigorous shaking or vortexing of the cell suspension.

**Note:** *iCell Astrocytes can be concentrated at thaw. Transfer the cell suspension to a 15 ml centrifuge tube and centrifuge at 300 x g for 5 minutes. Aspirate the supernatant, leaving 1 ml in the centrifuge tube, and resuspend the cell pellet in astrocyte culture medium.*

**Note:** *Thaw no more than 3 vials of iCell Astrocytes at one time. Once thawed, combine the contents of the cryovials before adding the rinse and final volume of astrocyte culture medium. Follow the timing outlined in steps 6 and 7.*

## Chapter 6. Plating iCell Astrocytes

Notes

The recommended plating density for iCell Astrocytes is  $\sim 4.0 - 5.5 \times 10^4$  viable cells/cm<sup>2</sup>.

1. Remove a sample of cells to perform a cell count using a hemocytometer (using trypan blue exclusion to identify viable cells) or an automated cell counter.
2. Dilute the cell suspension using room temperature astrocyte culture medium to obtain a desired cell plating density.
3. Aspirate the Matrigel from pre-coated plates and immediately dispense the cell suspension.
4. Culture iCell Astrocytes in a cell culture incubator at 37°C, 5% CO<sub>2</sub>.

### Expected Cell Density

$\sim 4.0 - 5.5 \times 10^4$  viable cells/cm<sup>2</sup> is the recommended starting density of iCell Astrocytes for most cell-based assays. However, the optimal density of iCell Astrocytes per unit of surface area can be assay dependent and must be determined empirically based on the intended use. The following table provides the desired cell number and plating volume for several common culture vessels.

| Culture Vessel             | Surface Area (cm <sup>2</sup> ) | Plating Volume (ml) | Cell Number ( $\sim 4.0 - 5.5 \times 10^4$ cells/cm <sup>2</sup> ) |
|----------------------------|---------------------------------|---------------------|--|
| 6-well Cell Culture Plate  | 9.6                             | 3                   | $38 - 52 \times 10^4$  |
| 24-well Cell Culture Plate | 1.9                             | 0.6                 | $7.6 - 10.4 \times 10^4$   |
| 96-well Cell Culture Plate | 0.32                            | 0.2                 | $1.2 - 1.7 \times 10^4$  |

**Table 2: Summary of Recommended Volumes and Measures**

*All volumes and measures are per well.*

## Chapter 7. Maintaining iCell Astrocytes

1. Immediately before use, equilibrate the astrocyte culture medium in a 37°C water bath.
2. 24 hours post-plating iCell Astrocytes, carefully aspirate the spent medium and replace (100% exchange) with the appropriate volume of astrocyte culture medium. Recommend volumes are as follows:
  - **6-well cell culture plate:** 2 ml/well
  - **24-well cell culture plate:** 0.6 ml/well
  - **96-well cell culture plate:** 200 µl/well
3. Replace the spent medium every 2 - 3 days.
4. Culture iCell Astrocytes in a cell culture incubator at 37°C, 5% CO<sub>2</sub>.

# Appendices

---

## Appendix A. Intellectual Property Rights, Use Restrictions, and Limited License

**A. OWNERSHIP.** The Products are covered by pending patents and patents: [www.cellulardynamics.com/patents](http://www.cellulardynamics.com/patents). Customer has a limited license to use the Products for internal research purposes for the sole benefit of the Customer, subject to the use restrictions and third party licenses included in subsections B and C of this Appendix A. Customer acknowledges and agrees that the receipt or purchase of the Products by Customer shall not be construed as a transfer of any title or the grant of any rights in or to the intellectual property embodied in the Products owned or licensed by Cellular Dynamics. In particular, no right or license to make, have made, offer to sell, or sell the Products, to modify or reproduce the Product or any part thereof, or to use the Products in combination with any other product(s), except product(s) provided or expressly licensed to Customer by Cellular Dynamics for such use, is implied or conveyed by the sale or transfer of Products to Customer.

**B. USE RESTRICTIONS.** The Products are licensed for internal research purposes only, and may not be used for any other purpose. The Products must be used in accordance with this User's Guide to which Customer, by ordering and accepting the Products, agrees. Customer shall not make, have made, offer to sell, or sell the Products. Customer shall not use the Products (or any modifications Customer makes to the Products or any cells derived, developed or expanded from the Products) in (i) the manufacture of any products, or (ii) any services for a third party. Customer may not transfer the Products (or any modifications Customer makes to the Products or any cells derived, developed or expanded from the Products) to any third party without Cellular Dynamics' prior written consent. Customer shall not reverse engineer the Products. Customer shall not use the Products, components or modifications thereof, or any cells derived, developed or expanded therefrom, in humans, in clinical trials, for diagnostic purposes involving human subjects, or for any investigational or other therapeutic use. Customer shall not use the Products directly or indirectly to derive or make any human gamete or gamete precursor cell. Customer shall use the Products in accordance with all applicable laws and regulations and any applicable institutional review board approved protocol and/or privacy office approval. Customer is not entitled to receive any data or information from Cellular Dynamics that directly identifies the donor of the biological materials from which the Products indirectly are derived or were made. Customer shall not attempt in any way to determine the identity of the donor of the biological materials from which the Products indirectly are derived or was made.

**C. DATA.** Customer agrees that if described on Customer's product quotation from Cellular Dynamics it will provide data and information as described therein to Cellular Dynamics regarding Customer's use of the Products.

---

## Appendix B. Product Provided "AS IS"

A. The Products are sold or provided "AS IS."

B. TO THE FULLEST EXTENT PERMITTED BY APPLICABLE LAW, CELLULAR DYNAMICS DISCLAIMS, ALL REPRESENTATIONS, AND WARRANTIES, EXPRESS OR IMPLIED (INCLUDING BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE OR NON-INFRINGEMENT), AND ALL LIABILITY, WITH RESPECT TO THE PRODUCTS; AND BY ORDERING AND ACCEPTING THE PRODUCTS CUSTOMER WAIVES ALL RIGHTS AND REMEDIES OTHERWISE AVAILABLE WITH RESPECT TO THE PRODUCTS.

C. Customer will be solely responsible for (i) Customer's use of the Products for a purpose or in a manner other than that for which they were designed or that is permitted or in breach of the Use Restrictions above; (ii) Customer's failure to follow this User's Guide for the use, storage, and handling of the Products however such failure is caused; (iii) Customer's failure to comply with any of the provisions of Appendix A above; and (iv) any abuse, other misuse



or neglect of the Products by Customer or any damage or loss of the Products by events or occurrences beyond a person's (e.g., Cellular Dynamics') control including without limitation, accident, fire, vandalism and natural disasters (acts of God).

D. Customer acknowledges and agrees that Cellular Dynamics may fill Customer's order with any number of units of Products. Such units may be more units than Customer ordered. Customer will not be charged extra for any adjustments made by Cellular Dynamics. Because the number of cells in a unit may vary from lot to lot, Cellular Dynamics reserves the right to fill the order with that number of units which is sufficient to fill Customer's order.

---

## **Appendix C. Limited Liability**

TO THE FULLEST EXTENT PERMITTED UNDER APPLICABLE LAW, CELLULAR DYNAMICS SHALL NOT HAVE ANY LIABILITY FOR INCIDENTAL, COMPENSATORY, PUNITIVE, CONSEQUENTIAL, INDIRECT, SPECIAL OR OTHER SIMILAR DAMAGES, HOWEVER CAUSED AND REGARDLESS OF FORM OF ACTION WHETHER IN CONTRACT, TORT (INCLUDING NEGLIGENCE), STRICT PRODUCT LIABILITY OR OTHERWISE, EVEN IF CELLULAR DYNAMICS HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. NOTWITHSTANDING ANY OTHER TERM OR IMPLICATION TO THE CONTRARY, UNDER NO CIRCUMSTANCES SHALL CELLULAR DYNAMICS' LIABILITY TO CUSTOMER EXCEED THE AMOUNT PAID BY CUSTOMER FOR THE PRODUCTS TO CELLULAR DYNAMICS.

© 2014 Cellular Dynamics International, Inc. All rights reserved.